



# TGEX-HC-mG2b-Zeo Expression Vector

## INSTRUCTION MANUAL

---

TGEX-HC-mG2b-Zeo Transient Mammalian Expression Vector

**Catalog #:** MX034

**Version:** A1.2 – February 2024



---

# Table of Contents

Limited Use License for TGEX Vector Series	4
<b>Description</b>	<b>5</b>
Introduction	5
Content, Shipping & Storage	5
Limited Product Warranty	5
TGEX Vector Series	6
Vector Map	6
Cloning Site	7
Feature Table	8
Restriction Site Summary	8
<b>Experimental Procedures</b>	<b>11</b>
General Molecular Biology Techniques	11
Plasmid Maintenance	11
Cloning into TGEX-HC-mG2b	11
Sequencing of Inserts	12
Antibody Expression	12
<b>Appendix</b>	<b>14</b>
MSDS Information	14
Quality Control	14
Technical Support	14
References	14

---

## Limited Use License for the TGEX Vector Series

---

As a condition of sale of this product to you, and prior to using this product, you must agree to the terms and conditions of this license. Antibody Design Labs grants to the buyer with the sale of any of its **TGEX™** vectors (the “Product”) a non-exclusive, non-transferable and limited license to use the Product in research only conducted by the buyer. Such license specifically excludes the right to sell or otherwise transfer the Product, its components or derivatives thereof to third parties. No modifications to the Product may be made without express written permission from Antibody Design Labs. The buyer is not granted a license to use the Product for human or animal therapeutic, diagnostic, or prophylactic purposes.

Antibody Design Labs does not warrant that the use or sale of the Product, the use thereof in combination with other products, or the use of the Product in the operation of any process will not infringe the claims of any United States or other patent(s).

If the buyer is not willing to accept the limitations of this license, without modification, buyer may refuse this license by returning the Product unopened and unused. By keeping or using the Product, buyer implicitly agrees to be bound by the terms of this license.

Entities wishing to use the Product for commercial purposes are required to obtain a license from Antibody Design Labs. Commercial purposes may include, but are not limited to: use of the Product in manufacturing, use of the Product to provide a service, use of the Product for therapeutic or diagnostic purposes, or resale of the Product, whether or not such Product is resold for use in research. For information on purchasing a commercial license to the Product, please contact a licensing representative by phone at (858) 480-6213 or by e-mail at [info@abdesignlabs.com](mailto:info@abdesignlabs.com).

All trademarks are the property of their respective owners.

# Description

## Introduction

The **TGEX™** vector series is designed for the rapid expression of antibody molecules by transient gene expression in mammalian cells in suspension culture. This new series, version 7, also allows the selection of stable transformants through the use of Zeocin®. This vector series features a cytomegalovirus (CMV) promoter, the adenovirus tripartite leader sequence (TPL) (Logan 1984, Mariati 2010), a composite leader sequence (hybrid between a mammalian kappa leader and the bacterial pelB leader (Valadon 2006)) followed by a universal double-inverted BsaI cloning site to insert antibody variable regions. The constant regions are derived from species-specific IgG, kappa and lambda sequences. The 3' end on the transcription unit is composed of an IRES element, the Zeocin-resistance gene (Sh ble from *Streptoalloteichus hindustanus*), a Woodchuck hepatitis virus post-transcriptional regulatory element, and the rabbit beta-globin polyadenylation signal. Transfection of the **TGEX™** expression vectors harboring antibody variable region inserts in widely available cell lines using large-scale transfection technologies (see experimental procedure) typically yields antibody titers between 50 and 250 mg/L in serum-free conditions in just a few days.

The **TGEX™-HC-mG2b-Zeo** vector is designed for the expression of a heavy chain variable region with the constant region of the murine IgG2b heavy chain. Expression of full-length antibody molecules is achieved by co-transfection with a light chain variable region cloned into one the **TGEX™** vectors, see below examples to achieve varied antibody formats.

### Combination of vectors to desired antibody format (examples)

COMBINATION	FORMAT	PURIFICATION
TGEX™-HC-hG1-Zeo + TGEX™-LC-hK-Zeo	Full length human or chimeric human IgG1/K	Protein A or G
TGEX™-FH-hG1-Zeo + TGEX™-LC-hK-Zeo	Human or chimeric IgG1/K Fab fragment	Protein L, G, or IMAC
TGEX™-FC-hG1-Zeo	Human IgG1 Fc fusion	Protein A or G
TGEX™-SCblue-Zeo	Human IgG1 scFv-Fc fusion	Protein A or G

## Content, Shipping & Storage

### Content

VECTOR	COMPOSITION	AMOUNT
TGEX-HC-mG2b-Zeo	20 µl at 0.5 µg/µl of DNA vector in DNA Conservation Buffer (Tris-HCL 5 mM, EDTA 0.1 mM, pH 8.5)	10 µg

### Shipping & Storage

**TGEX-HC-mG2b-Zeo** vector is shipped on wet ice. Upon receipt, store the vector at -20°C.

## Limited Product Warranty

This warranty limits our liability to the replacement of this product. No other warranties of any kind express or implied, including, without limitation, implied warranties of merchantability or fitness for a particular purpose, are provided by Antibody Design Labs. Antibody Design Labs shall have no liability for any direct, indirect, consequential, or incidental damages arising out of the use, the results of use, or the inability to use this product.

**For research use only; not intended for any animal or human therapeutic or diagnostic use.**

## TGEX™ Vector Series

### TGEX™ vector series for desired isotypes

VECTOR	CATALOG	ISOTYPE	DESCRIPTION
TGEX™-HC-hG1-Zeo	MX026	Human IgG1	Heavy chain expression plasmid for human IgG1
TGEX™-FH-hG1-Zeo	MX023	Human IgG1 CH1	For the expression of human IgG1 Fab fragments
TGEX™-FC-hG1-Zeo	MX025	Human IgG1 Fc	For the expression of human IgG1 Fc fusions
TGEX™-HC-hG2-Zeo	MX027	Human IgG2	Heavy chain expression plasmid for human IgG2
TGEX™-HC-hG3-Zeo	MX028	Human IgG3	Heavy chain expression plasmid for human IgG3
TGEX™-HC-hG4[S228P]-Zeo	MX029	Human IgG4[S228P]	Heavy chain expression plasmid for human IgG4
TGEX™-LC-hK-Zeo	MX030	Human Kappa	Light chain expression plasmid for human Kappa
TGEX™-LC-hL2-Zeo	MX031	Human Lambda 2	Light chain expression plasmid for human Lambda 2
TGEX™-HC-mG1-Zeo	MX032	Murine IgG1	Heavy chain expression plasmid for murine IgG1
TGEX™-HC-mG2a-Zeo	MX033	Murine IgG2a	Heavy chain expression plasmid for murine IgG2a
TGEX™-HC-mG2b-Zeo	MX034	Murine IgG2b	Heavy chain expression plasmid for murine IgG2b
TGEX™-HC-mG3-Zeo	MX035	Murine IgG3	Heavy chain expression plasmid for murine IgG3
TGEX™-LC-mK-Zeo	MX036	Murine Kappa	Light chain expression plasmid for murine Kappa
TGEX™-LC-mL1-Zeo	MX037	Murine Lambda 1	Light chain expression plasmid for murine Lambda 1
TGEX™-LC-mL2-Zeo	MX038	Murine Lambda 2	Light chain expression plasmid for murine Lambda 2
TGEX™-HC-rbG-Zeo	MX039	Rabbit IgG	Heavy chain expression plasmid for rabbit IgG
TGEX™-LC-rbKb4-Zeo	MX040	Rabbit Kappa	Light chain expression plasmid for rabbit Kappa
TGEX™-HC-dG1-Zeo	MX046	Dog IgG1	Heavy chain expression plasmid for dog IgG1
TGEX™-HC-dG2-Zeo	MX047	Dog IgG2	Heavy chain expression plasmid for dog IgG2
TGEX™-HC-dG3-Zeo	MX048	Dog IgG3	Heavy chain expression plasmid for dog IgG3
TGEX™-HC-dG4-Zeo	MX049	Dog IgG4	Heavy chain expression plasmid for dog IgG4
TGEX™-LC-dK-Zeo	MX050	Dog Kappa	Light chain expression plasmid for dog Kappa
TGEX™-LC-dL-Zeo	MX051	Dog Lambda	Light chain expression plasmid for dog Lambda

### TGEX™ vector series for Fc-engineered antibodies


VECTOR	CATALOG	ISOTYPE	DESCRIPTION
TGEX™-HC-hG1[EA]-Zeo	MX041	Human IgG1	Human IgG1 heavy chain with increased ADCC/CDC <i>in vitro</i>
TGEX™-HC-hG1[NA]-Zeo	MX042	Human IgG1	Human IgG1 aglycosylated heavy chain
TGEX™-HC-hG1[LALA-PG]-Zeo	MX043	Human IgG1	Human IgG1 heavy chain with decreased ADCC/CDC <i>in vitro</i>
TGEX™-HC-hG1[YTE-KF]-Zeo	MX044	Human IgG1	Human IgG1 heavy chain with increased serum half-life
TGEX™-HC-hG4[SPLP-PG]-Zeo	MX045	Human IgG4[S228P]	Human IgG4 heavy chain with decreased ADCC/CDC <i>in vitro</i>

### TGEX™ control vector, universal expression and Fc fusions

VECTOR	CATALOG	USE	DESCRIPTION
TGEX™-AC-Zeo	MX020	Any expressions	Universal expression vector
TGEX™-eGFP-Zeo	MX022	Transfection	Control plasmid for monitoring transient transfections
TGEX™-SCblue-Zeo	MX024	scFv cloning vector	For the transfer of scFv from any PADL phagemid vector and expression as an scFv-Fc fusion

## Vector Map

The figure below illustrates the main features of **TGEX-HC-mG2b-Zeo** expression vector. The full vector sequence is available online for download in varied formats on the product web page; the total length of the vector is 5846 bp.

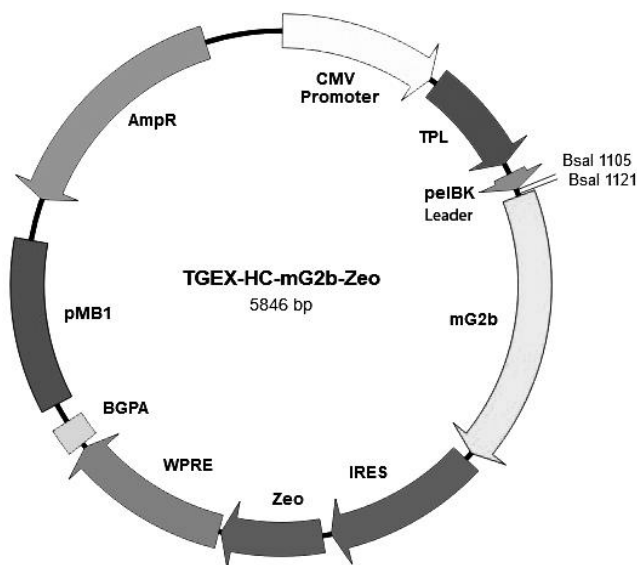

## Cloning Site

Following is an illustration of **TGEX-HC-mG2b-Zeo** cloning site from the EcoRI site and onward. The VH domain is inserted in a double inverted Bsal cloning site located between the pelBK leader sequence and the constant regions. The four base pair overhangs after a restriction digestion with Bsal are boxed and grayed.

```

          tgex-S3                                     EcoRI
          |----->                                     |
971  GAAAGGCGTC TAACCAGTCA CAGTCGCAAG TTTAAACGGA TCTCTAGCGA
          |
          | pelBK leader sequence
          | MetGluThrA spThrLeuLe uLeuTrpVal
1021  ATTCGGCTTG GGCCGCCACC ATGGAGACAG ACACACTCCT GCTATGGGTA
          |
          | Sfil                                     <Bsal                                     MluI
          | |                                     | |                                     |
          | LeuLeuLeuL euAlaAlaGl nProAlaMet Ala
1071  CTGCTGCTCT TAGCGGCCCA GCCGGCCATG GCA GGAGACC TAACGCGTAT
          |
          | >Bsal                                     mG2b
          | |----->
          | Ala LysThrThrP roProSerVa lTyrProLeu AlaProG...
1121  GGTCTCAGCT AAAACAACAC CCCCATCAGT CTATCCACTG GCCCCTGGGT
  
```

## Feature Table

The features of **TGEX-HC-mG2b-Zeo** transient expression vector are highlighted in the following table.

FEATURE	LOCATION	DESCRIPTION
Promoter	5-585	CMV promoter.
TPL	612-1000	Adenovirus tripartite leader sequence (Logan 1984, Mariati 2010).
pelBK leader	1041-1103	Hybrid kappa/pelB leader peptide sequence. The cleavage occurs on the C-terminal side of the terminal alanine.
Murine IgG2b CDS	1128-2135	Sequence encoding the murine IgG2b heavy chain sequence comprising the CH1, CH2, and CH3 domains with the hinge region. The sequence is intronless; the C-terminal lysine residue is encoded.
IRES	2164-2737	Internal Ribosome Entry Site.
Zeo	2771-3145	Sh ble gene from <i>Streptoalloteichus hindustanus</i> conferring resistance to Zeocin.
WPRE	3154-3742	Woodchuck hepatitis virus post-transcriptional regulatory element.
BGpA	3758-3856	Rabbit beta-globin polyadenylation signal sequence.
pMB1 origin	3930-4549	pBR322 origin for replication in <i>E. coli</i> with a temperature-sensitive high copy-number phenotype (Lin-Chao 1992).
TEM1 beta-lactamase	5564-4704	Ampicillin resistance for selection in <i>E. coli</i> .

## Restriction Site Summary

Enzyme	Site	Nb	Position	Strand	Isoschizomers
AbsI	CC <sup>^</sup> TCGAGG	1	1493		
AlfI	(10/12) GCANNNNNNTGC (12/10)	1	3108		
AloI	(7/12) GAACNNNNNTCC (12/7)	1	1250		
ArsI	(8/13) GACNNNNNTTYG (11/6)	1	851		
AvrII	C <sup>^</sup> CTAGG	1	2315		AspA2I BlnI XmaJI
BcgI	(10/12) CGANNNNNNTGC (12/10)	1	5293		
BsaXI	(9/12) ACNNNNNTCC (10/7)	1	3007	-	
BsePI	G <sup>^</sup> CGCGC	1	2807		BssHII PauI PteI
Bsp1407I	T <sup>^</sup> GTACA	1	5808		BsrGI BstAUI
BstEII	G <sup>^</sup> GTNACC	1	2122		BstPI Eco91I EcoO65I PspEI
BstZ17I	GTA <sup>^</sup> TAC	1	1836		BssNAI Bst1107I
CspCI	(11/13) CAANNNNGTGG (12/10)	1	407		
EcoNI	CCTNN <sup>^</sup> NNNAGG	1	1208		BstENI XagI
EcoRI	G <sup>^</sup> AATTC	1	1019		
FseI	GGCCGG <sup>^</sup> CC	1	3042		RigI
FspI	TGC <sup>^</sup> GCA	1	4997		Acc16I NsbI
KpnI	GGTAC <sup>^</sup> C	1	2605		Acc65I Asp718I
MauBI	CG <sup>^</sup> CGCGCG	1	2806		
MluI	A <sup>^</sup> CGCGT	1	1113		
NheI	G <sup>^</sup> CTAGC	1	2142		AsuNHI BmtI BspOI
NotI	GC <sup>^</sup> GGCCGC	1	3146		CciNI
PasI	CC <sup>^</sup> CWGGG	1	1163		
PciI	A <sup>^</sup> CATGT	1	2654		PscI
PmeI	GTTT <sup>^</sup> AAAC	1	1000		MssI
PspXI	VC <sup>^</sup> TCGAGB	1	1493		
PvuI	CGAT <sup>^</sup> CG	1	5144		Ple19I
SalI	G <sup>^</sup> TCGAC	1	2765		
SexAI	A <sup>^</sup> CCWGGT	1	2934		CsiI MabI



SfiI	GGCCNNNN^NGGCC	1	1085	
SgrAI	CR^CCGGYG	1	2884	
SnaBI	TAC^GTA	1	357	BstSNI Eco105I
SpeI	A^CTAGT	1	18	AhlI BcuI
XbaI	T^CTAGA	1	2759	
AarI	CACCTGC (4/8)	2	1367	PaqCI
		2	2493	
ApaI	GGGCC^C	2	1428	Bsp120I PspOMI
		2	2277	
BamHI	G^GATCC	2	1257	
		2	2753	
BglII	A^GATCT	2	3873	
		2	5697	
Bpu10I	CCTNAGC (-5/-2)	2	930	
		2	1227	
BsmI	GAATGC (1/-1)	2	2310	- Mva1269I PctI
		2	2343	
BspMI	ACCTGC (4/8)	2	1368	Acc36I BfuAI BveI
		2	2494	
BtrI	CACGTC (-3/-3)	2	2705	- AjiI BmgBI
		2	2903	
DraIII	CACNNN^GTG	2	2522	AdeI
		2	3119	
EagI	C^GGCCG	2	3102	BseX3I BstZI EclXI Eco52I
		2	3147	
Eam1105I	GACNNN^NNGTC	2	1504	AhdI BmeRI DriI
		2	4772	
Eco31I	GGTCTC (1/5)	2	1105	- Bso31I BsaI BspTNI
		2	1121	
FalI	(8/13)AAGNNNNCTT (13/8)	2	784	
		2	1215	
HindIII	A^AGCTT	2	2736	
		2	3879	
NdeI	CA^TATG	2	252	FauNDI
		2	3835	
OliI	CACNN^NNGTG	2	1397	AleI
		2	1571	
PmaCI	CAC^GTG	2	1571	AcvI BbrPI Eco72I PmlI PspCI
		2	2478	
PstI	CTGCA^G	2	1296	BspMAI
		2	1370	
PvuII	CAG^CTG	2	642	
		2	1619	
SacI	GAGCT^C	2	583	Ecl136II EcoICRI Eco53kI
				Psp124BI SstI
		2	1825	
SacII	CCGC^GG	2	740	Sfr303I KspI SgrBI Cfr42I
		2	3654	
SmaI	CCC^GGG	2	2749	Cfr9I TspMI XmaI
		2	2859	
VspI	AT^TAAT	2	25	AseI PshBI
		2	4948	
XhoI	C^TCGAG	2	966	Sfr274I Paer7I SlaI
		2	1494	
XmnI	GAANN^NNTTC	2	2372	Asp700I MroXI PdmI
		2	5372	

Absent Sites:

AanI, AccIII, AfeI, AflIII, AgeI, AjuI, Aor13HI, Aor51HI, AscI, AsiGI, AsiSI, AsuII, AxyI, BaeI, BarI, BbvCI, BclI, BfrI, BlpI, BoxI, BplI, Bpu1102I, Bpu14I, Bsa29I, BsaBI, Bse21I, Bse8I, BseAI, BseCI, BseJI, BsgI, BshTI, BshVI, BsiWI, BsmBI, Bsp119I, Bsp13I, Bsp1720I, Bsp68I, BspDI, BspEI, BspQI, BspT104I, BspTI, BstAFI, BstAPI, BstBI, BstPAI, BstXI, Bsu15I, Bsu36I, BsuTUI, BtuMI, ClaI, CpoI, CspAI, CspI, DinI, Eco147I, Eco32I, Eco47III, Eco81I, EcoRV, EcoT22I, EgeI,

---

EheI, Esp3I, FbaI, FspAI, HpaI, I-CeuI, I-PpoI, I-SceI, KasI, KflI, Kpn2I, Ksp22I, KspAI, LguI, MfeI, Mly113I, Mph1103I, MreI, MroI, MspCI, MunI, NarI, NruI, NsiI, NspV, PI-PspI, PI-SceI, PacI, PaeI, PalAI, PceI, PciSI, Pfl23II, PflFI, PinAI, PluTI, PshAI, PsiI, PspLI, PsrI, PsyI, RgaI, RruI, Rsr2I, RsrII, SapI, SbfI, SdaI, SfaAI, SfoI, SfuI, SgfI, SgrDI, SgsI, SmiI, SphI, SrfI, Sse8387I, SseBI, SspDI, StuI, SwaI, Tth111I, Vha464I, XcmI, Zsp2I.

---

# Experimental Procedures

---

## General Molecular Biology Techniques

---

Molecular biology should be conducted under the supervision of a qualified instructor trained to standard safety practice in a molecular biology laboratory environment. Standard molecular biology procedures can be found in a general molecular biology handbook such as Sambrook (1989).

---

## Plasmid Maintenance

---

Propagation and maintenance of TEGX vectors is obtained on any *recA1, endA1 E. coli* strain using LB or 2xYT medium supplemented with ampicillin (100 µg/ml) as a selection marker and incubated at 37°C with agitation. TEGX vectors are derivatives of pBR322 with a high copy number origin of replication and usually gives high yields of plasmid DNA with most standard laboratory strains such as XL1-blue or DH5α. The high copy number phenotype is temperature-sensitive and requires incubation at 37°C (Lin-Chao 1992). Some DNA stabilizing strains are known to produce smaller amounts of plasmid DNA. In case of issues, we recommend using XL10-Gold® from Agilent Technologies, Inc., on which TGEX plasmid DNA can be isolated in large quantities.

---

## Cloning into TGEX-HC-mG2b-Zeo

---

### *In Silico* Design

A double inverted Bsal cloning site separates the pelBK leader peptide sequence from the beginning of the antibody constant region. After digestion with Bsal, the heavy chain variable domain is inserted in-frame between the two elements as illustrated below. During the cloning, the two Bsal sites are eliminated.

```

                SfiI
                |
    LeuLeuLeuL euAlaAlaGl nProAlaMet Ala
1071 CTGCTGCTCT TAGCGGCCCA GCCGGCCATGGCA [VH] GCTAAAACAACAC CCCCATCAGT
                mG2b
                |
    Ala LysThrThrP roProSerVa

```

### Vector Digestion

Bsal alone is sufficient to open the cloning site. Please, consult the documentation of your restriction enzyme provider for optimal conditions.

### Cloning with FAST-Licase

Homologous recombination is the easiest and most efficient method to clone seamlessly antibodies into TGEX vectors. We recommend the **FAST-Licase™** (Antibody Design Labs cat# MB101S & MB101L) with the following overhangs. For synthetic dsDNA constructs, the overhangs are added on each side. For PCR-amplified constructs, primers should contain the overhang followed by the antibody priming area. The **FAST-Licase™** reaction contains the insert plus the purified vector digested with Bsal (see kit instructions).

pelBK overhang            5' – CGGCCCAGCCGGCCATGGCA

mG2b overhang            5' – GCTAAAACAACACCCCATC (reverse complement for PCR primers)

---

## Primer Design for Restriction Cloning with the BsaI Sites

Oligo1 is an example of primer designed to amplify a VH domain sequence and clone it into the BsaI site situated next to the pelBK leader. A minimum of 2 nucleotides is recommended to cut BsaI site close to the end (source New England Biolabs); these 2 nucleotides are followed by the BsaI site GGTCTC and the last five nucleotides of the leader sequence; the resulting NNGGTCTCTGGCA extension is added 5' to the VH domain primer.

**Oligo1**      5' -NNGGTCTCTGGCA-VH-start

Oligo2 is an example of primer designed to amplify the VH domain from the end of the J region including a BsaI site compatible with **TGEX-HC-mG2b-Zeo** cloning site.

BsaI is a type IIS restriction enzyme that cuts outside of its recognition site. The second BsaI site of the cloning site will be cut immediately before the start of the murine IgG2b constant region sequence, exactly 1 base after the end of the site and 5 bases further on the opposite strand, thus freeing a 5' 4-base overhang GCAT on the sense strand. In Oligo2, a BsaI recognition site is situated symmetrically to the vector and will generate a complementary overhang.

**Oligo2**      5' -NNGGTCTCGTAGC-JH-end

After digestion and ligation, all the two BsaI sites, from both insert and vector, will be removed, resulting in a scarless insertion of the antibody domain.

## Alternative to BsaI Sites

All restriction enzymes that generate 4-base long 5' overhangs can be used in place of BsaI; this is the case for example of BsmBI (CGTCTC(1/5)), another type IIS restriction enzyme. This option could come in handy when the VH domain contains another BsaI preventing cloning.

---

## Sequencing of Inserts

The following primers give a strong PCR amplification of the TGEX vector series inserts and the antibody constant regions between the EcoRI site and the NheI site. The primer tgex-S3 can be used to sequence the VH domain in full.

**tgex-S3**      5'- AGGCGTCTAACCAGTCACAGTC

**ires-R**      5'- GAATAAGGCCGGTGTGCGTT

---

## Antibody Expression

### Cell Lines

Cell lines adapted for culture in suspension and serum-free conditions are recommended. HEK293 and CHO cells are often used for antibody expression by transient transfection; you can either adapt your own cell line or obtain it from a supplier (e.g. Life Technology). HEK293 cells are particularly well suited for expression using **TGEX™** vector series.

### Transient Transfection

Many transfection reagents especially designed for transient transfection are commercially available from different providers (e.g. Life Technologies, Mirus Bio LLC). We recommend testing the transfection conditions with a reporter plasmid first to determine the percentage of cells effectively transfected and optimal transfection conditions; fluorescent reporters are often used with that purpose, e.g. TGEX™-eGFP-Zeo (Antibody Design Labs cat# MX022). Similarly, any condition known to boost expression should be carefully tested in your system before being scaled up. We did observe an increase in

---

expression in HEK293 cells upon exposure to sodium valproate (Backliwal 2008). Boosters and enhancers are often included in commercially available transfection kits.

## Stable Cell Line Selection

Zeocin can be used very effectively to select stable cell lines secreting antibodies in a few weeks. The following protocol are suggestions that require adjustment to your particular constructs:

### WEEK 1.

- Day 0: Transfect in duplicate cells in a 6-well plate containing 2 ml culture per well;
- Day 3-4: Expand each well in two wells with 2 ml culture per well.
- Day 6: Increase the volume to 4 ml per well and add Zeocin at 100 µg/ml.

### WEEK 2-3

- Maintain a good cell density and the Zeocin concentration while replenishing the culture with fresh medium as needed.

### WEEK 3-4

- Continue selection with possibly a higher Zeocin concentration (up to 1000 µg/ml).
- Proceed to single cloning and analyze stable transformants for expression.

## Light chain to Heavy chain ratio

We recommend starting with a 1:1 light chain to heavy chain ratio during transfection. We observed many antibodies with a better expression at a 2:1 light chain to heavy chain ratio although each antibody requires fine tuning for optimal expression.

---

# Appendix

---

## MSDS Information

---

Material Safety Data Sheets are available on Antibody Design Labs website at the corresponding product page.

---

## Quality Control

---

Specifications and quality control are detailed on the online product page. Antibody Design Labs certifies that the product will perform according to these specifications.

---

## Technical Support

---

Visit Antibody Design Labs website at [www.abdesignlabs.com](http://www.abdesignlabs.com) for technical resources, including manuals, vector maps and sequences, application notes, FAQs, etc.

For more information or technical assistance, call, write, or email us at:

Antibody Design Labs

4901 Morena Blvd, Suite 203

San Diego, CA 92117

Email: [support@abdesignlabs.com](mailto:support@abdesignlabs.com)

Phone: 1-877-223-3104 (Toll Free)

(Monday – Friday 9:00 AM – 5:00 PM PST)

---

## References

---

1. Logan J, & Shenk T. (1984). Adenovirus tripartite leader sequence enhances translation of mRNAs late after infection. *Proc Natl Acad Sci USA*, 81(12):3655–9.
2. Mariati, Ho SCL, Yap MGS, & Yang Y. (2010). Evaluating post-transcriptional regulatory elements for enhancing transient gene expression levels in CHO K1 and HEK293 cells. *Protein Expr Purif*, 69(1):9–15.
3. Valadon P, Garnett JD, Testa JE, Bauerle M, Oh P, & Schnitzer JE. (2006). Screening phage display libraries for organ-specific vascular immunotargeting in vivo. *Proc Natl Acad Sci USA*, 103(2):407–12.
4. Lin-Chao S, Chen WT, Wong TT (1992). High copy number of the PUC plasmid results from a ROM/ROP-suppressible point mutation in RNA II. *Mol Microbiol*, 6(22):3385–93.
5. Sambrook J, Fritsch EF, & Maniatis T. (1989). In *Molecular cloning: A laboratory manual*. Cold Spring Harbor Laboratory Press, NY, VOL. 1, 2, 3.
6. Backliwal G, Hildinger M, Kuettel I, Delegrange F, Hacker DI, Wurm FM. (2008). Valproic acid: A viable alternative to sodium butyrate for enhancing protein expression in mammalian cell cultures. *Biotechnol Bioeng*, 101(1):182–9.

*This product is subject to Antibody Design Labs Terms & Conditions of Sales available online at <http://www.abdesignlabs.com/terms/>.*

*© 2024 Antibody Design Labs. All rights reserved.*